

Tutorial

On this tutorial, we are going to demonstrate some features of the parKVFinder, including the graphical user interface (*parKVFinder PyMOL Tools*) and the command-line interface.

All files used on this tutorial can be found under **input** directory, on the **parKVFinder** directory.

parKVFinder PyMOL Tools

First, load **input/1FMO.pse** into PyMOL viewer, which loads two objects in your scene. The **1FMO** is a subunit of a protein kinase A and the **ligs_1FMO** is an adenosine (ADN) and a peptide kinase inhibitor (PKI).



Whole protein detection

The default parameters are designed to make a simple and fast whole protein detection.

On PyMOL, open **parKVFinder PyMOL Tools** under **Plugin** tab. The objects on the scene will be listed on the **Input PDB** listbox, on the **Main** tab. If not, press the **Refresh List**

The **Input PDB** selection sets which object will be analyzed by parKVFinder. Select **1FMO** on the listbox.



To run parKVFinder with the default parameters, just click **Run parKVFinder** button or press **Enter**.



After execution is complete, cavities PDB is loaded into PyMOL viewer as <Output Base Name>.KVFinder.output object and the results file is loaded on the **Results Visualization** tab. In addition, the focus automatically shifts to **Results Visualization** tab.



You can select cavities in the **Volume** or **Surface Area** lists to highlight them on a new object called **cavities**, helping to identify each cavity. Also, you can select cavities in the **Interface Residues**

list to highlight the residues around the cavities on a new object named **residues**.



Changing cavity ceiling

parKVFinder is all about parameter customization. One of parKVFinder's most powerful assets is the ability to manually set the cavity ceiling. parKVFinder works with a double probe system. A smaller probe, called Probe In, and a bigger one, called Probe Out, that defines two molecular surfaces with different molecular accessibility. The space left between these surfaces is considered cavities.

Let's show the effect of varying **Probe Out** and **Removal Distance** on the cavity ceiling.

First, we should copy the adenosine to a new object using the following PyMOL commands:

```
# Copy adenosine
select resn ADN
create adenosine, sele
delete sele
```

Also, copy the adenosine cavity (KAF) to a new object (adnsite) to compare the cavity ceiling from the previous execution.

```
# Copy adenosine cavity
select resn KAF and output.KVFinder.output
create adnsite, sele
delete sele
```

Finally, prepare the new PyMOL scene.

```
# Prepare PyMOL scene
color magenta, adnsite
disable
enable (adnsite, adenosine)
```



Adjusting Probe Out

As mentioned above, adjusting the Probe Out size changes the level of the cavity ceiling. So let's go back on the **Main** tab and change the **Probe Out** size to 8.0 Å. Run parKVFinder again.



Again, copy the adenosine cavity (KAF) to a new object (adnsitePO).

```

# Copy new adenosine cavity
select resn KAF and output.KVFinder.output
create adnsitePO, sele
delete sele

```

Finally, prepare the PyMOL scene.

```

# Prepare PyMOL scene
disable
enable (adenosine, adnsite, adnsitePO)


```



Note that the adenosine cavity detected with the 4 Å Probe Out (magenta) has a lower ceiling than that detected with the 8 Å probe (white). Therefore, by increasing the size of the Probe Out, the cavity ceiling is also raised.

Adjusting Removal Distance

Besides adjusting the Probe Out size, you can also adjust the Removal Distance to change the cavity ceiling. So let's go back to the **Main** tab and change the **Removal Distance** to 1.2 Å. Run parKVFinder again.

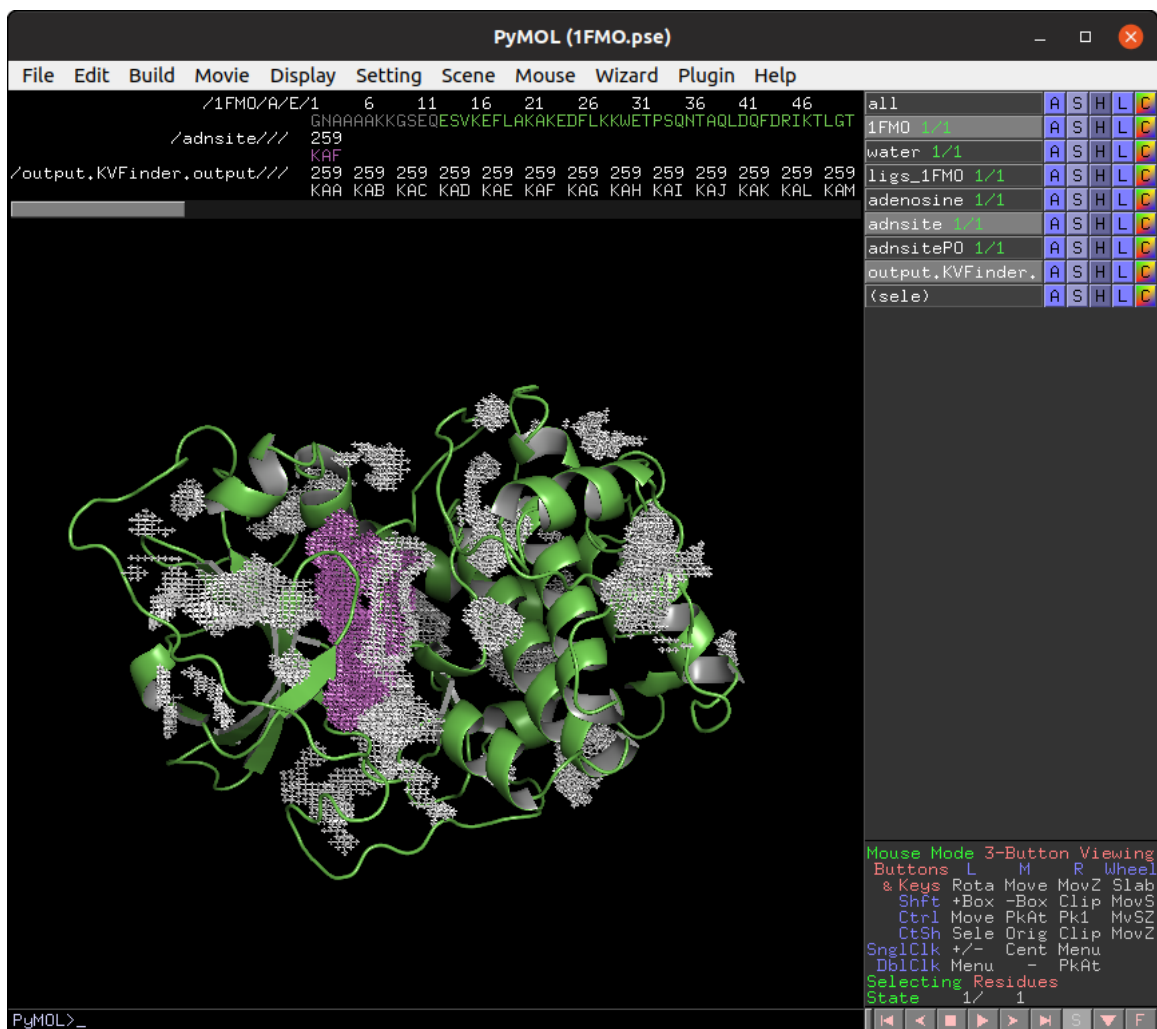


Again, copy the adenosine cavity (KAH) to a new object (adnsiteRD).

```
# Copy new adenosine cavity
select resn KAH and output.KVFinder.output
create adnsiteRD, sele
delete sele
```

Finally, prepare the PyMOL scene.

```
# Prepare PyMOL scene
disable
enable (adenosine, adnsite, adnsiteRD)
```

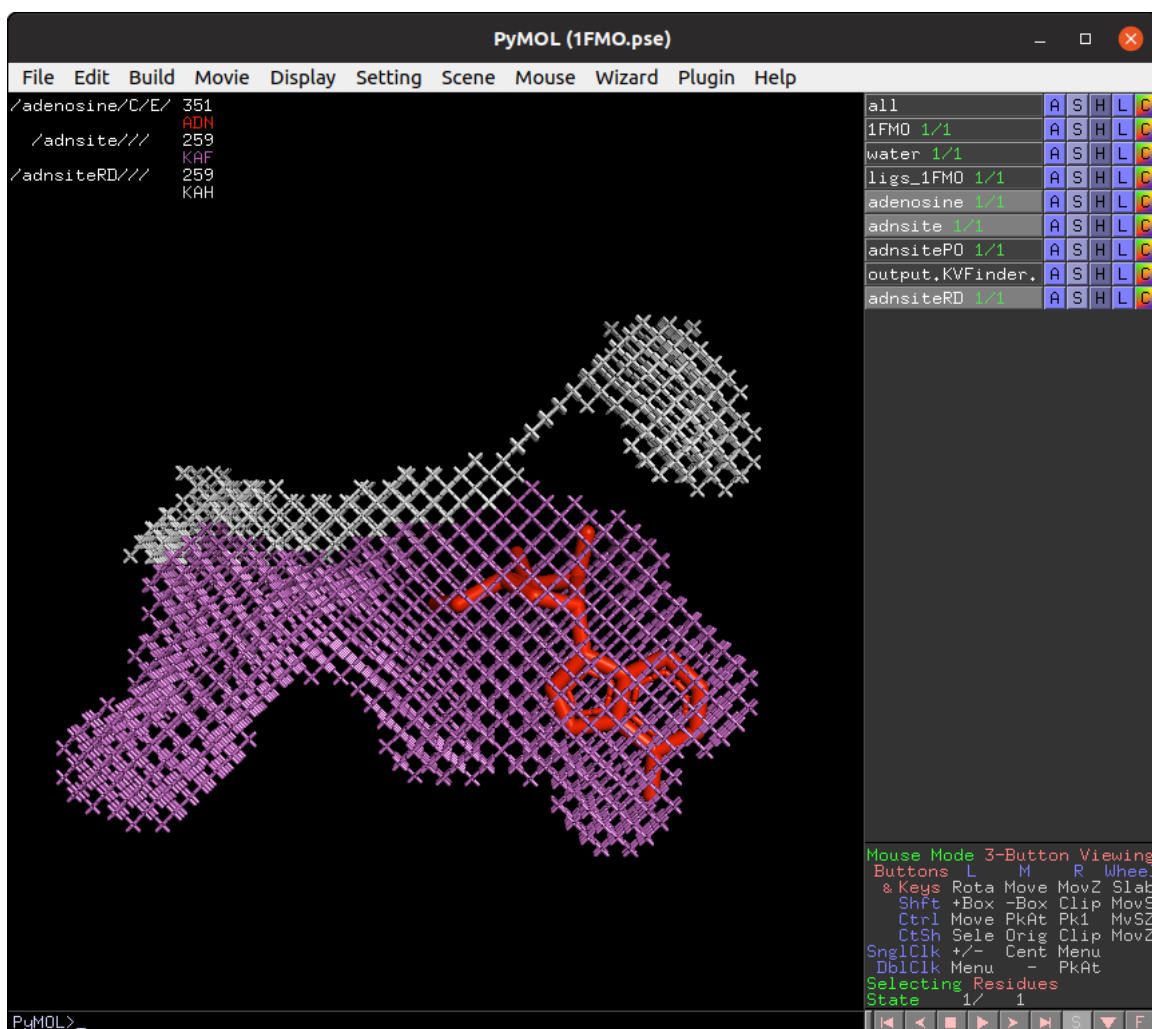Note that the adenosine cavity detected with the 2.4 Å Removal Distance (magenta) has a lower ceiling than that detected with the 1.2 Å (white). Therefore, by decreasing the Removal Distance, the cavity ceiling is also raised.

Furthermore, changing the cavity ceiling by varying Probe Out and Removal Distance also affects cavity segregation.

Note: Usually the Removal Distance adjustment is less time consuming than the Probe Out adjustment for similar effects.

Exploring enclosed regions with box adjustment mode

On Box Adjustment Mode, we can explore the sterred

On **Search Space** tab,

Cavity Segmentation - One of the main features of KVFinder is the ability to segment cavities in sub-pockets. In this tutorial we are going to demonstrate how to do it. For a meaningful biological example, please refer the the KVFinder paper.

Ligand adjustment mode

Command line interface

parKVFinder has a command-line interface, which can be useful for molecular dynamics and high-throughput analysis. It also handles the same parameters available in parKVFinder PyMOL Tools, except for box rotations in box adjustment mode.
